ordering information

Unique portfolio offering low & standard densities as well as high, standard & low viscosities

Squidperi

Reference	Viscosity	Density	Descripton
SQUIDPERI 34	34 cP	/	High viscosity & standard radiopacity
SQUIDPERI 34LD	34 cP	30% less than Squidperi 34	High viscosity & less radiopacity
SQUIDPERI 18	18 cP	/	Standard viscosity & standard radiopacity
SQUIDPERI 18LD	18 cP	30% less than Squidperi 18	Standard viscosity & less radiopacity
SQUIDPERI 12	12 cP	/	Low viscosity & standard radiopacity
SQUIDPERI 12LD	12 cP	30% less than Squidperi 12	Low viscosity & less radiopacity

squidperi compatibilities*



SonicDetachable tip microcathete



Eclipse 2L

Double lumen balloon microcathete



Copernic 2L

Double lumen balloon microcathete

*all DMSO compatible microcatheter

SQUIDPERI is intended to the embolisation of lesions in the peripheral vasculature, including arteriovenous malformations and hypervascular tumors. Class Ilb CE0297 in compliance with Medical Device Directive (MDD 93/42/EEC amended by 2007/47/EC). Manufactured by EMBO-FLÜSSIGKEITEN A.G. Carefully read the instructions for use before use. First CE marking: 2013. The content of this document, in particular data, information, trademarks and logos are BALT SAS's sole property © 2018 BALT and its affiliates, all rights reserved. All representation and/or reproduction, whether in part or in full, is forbidden and would be considered a violation of BALT SAS and its affiliates' copyrights and other intellectual proprietary rights. Any complete or partial reproduction, representation, modification, publication, transmission or distortion of the services or the contents of this document, whoever owns it, is forbidden, to the exception of a single copy on a computer reserved to the exclusive use of the copyist. This document with associated pictures are non-contractual and are solely dedicated to healthcare professionals and BALT's distributors (BALT's supplier's distributors). The products commercialized by BALT SAS shall exclusively be used in accordance with the instructions for use included in the boxes. DC033GB (03/2018)

Balt S.A.S

10, rue de la Croix Vigneron, 95160 Montmorency France Tél. : +33 (0)1 39 89 46 41 Tax : +33 (0)1 34 17 03 46

Fax: +33 (0)1 34 www.halt.fr



squidperi

the power of visibility, the choice of fluidness





squidperi

Last generation of EVOH co-polymer for peripheral usage

Squidperi is indicated in the embolisation of lesions in the peripheral vasculature, including arteriovenous malformations and hypervascular tumours.

Viscosities of

12 cP, 18 cP & 34 cP available

 The right fluidness & thickness, for the adapted procedure

12 cP exclusive formula penetration into the nidus &

enable deeper penetration into the nidus & to reach distal microvessels

18 cP regular formula


34 cP newcome For high flow pathologies Squidperi 34 Squidperi 34LD









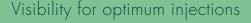

Homogenous radiopacity

LD UNIQUE formulas


To assess the angioarchitecture of the AVM (30% less TA grains) & the volume of liquid embolic injected



Micronization process

To minimize Tantalum powder grain siz

Cohesive penetration

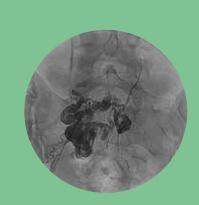

15 minutes after shaking

Squidperi

Micronized grain size of Tantalum powder

Precipitation of — Tantalum powder twice slower

Other liquid embolic agent

Siandara grain size o Tantalum powde

Regular
precipitation of
Tantalum powder

Each Squidperi kit contains

One 1,5 mL vial of Squidperi

One 1.5 mL vial of DMSC

One 1 cc Blue syringe for DMSC

Two 1 cc White syringes for Squidper

Two syringe adapters

