

Laryngeal and Esophageal Products

Har-El Pharyngeal Tube

The Har-El Pharyngeal Salivary Bypass Tube allows the surgeon to handle a variety of indications for treatment of esophageal fistulae. The proximal end features an indented funnel shape designed for anchoring the device at the level of the tongue base.

The tube has applications for: 1) the treatment of orocutaneous and pharyngocutaneous fistulae following head and neck oncologic resection, 2) placement in patients with an existing fistula prior to reconstructive procedures, and 3) prophylactic placement in patients who are at high risk for development of a fistula.

Intended Use

To facilitate the management of esophageal fistulae resulting from a variety of causes including surgery, radiation, advanced malignancy, trauma, and ingestion of caustic agents

Supplied


Sterile

For single use

1 Tube, 1 IFU per carton

Features

- Implant-grade silicone
- Non-adherent smooth surface
- Available in clear or radiopaque material

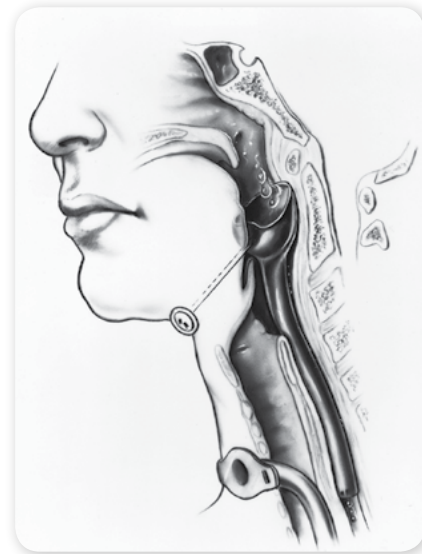

REF	
Clear	Radiopaque
355120	355120R

